INSTRUCTIONS FOR USE

Ambu[®] AuraFlex[™]

Single Use Laryngeal Mask - Sterile.


For use by medical professionals trained in airway management only.


Abridged version


Ambu





Symbol Indication	MD	MR	STERILE R	®	∕ÇN .	UK ÇA	UK RP	
en	Medical Device	MR unsafe	Sterilized using irradiation Single sterile barrier system	Do not use if the product sterilization barrier or its packaging is damaged	Country of manufacturer	UK Conformity Assessed	UK Responsible Person	Importer (For products imported into Great Britain only)

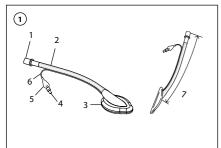

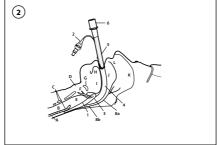

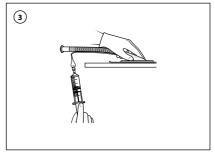

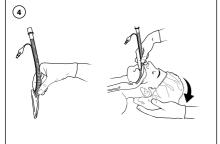

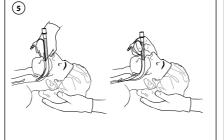



Content Page Content

English.......10-14

Page

1.1. Intended use/Indications for use

The Ambu AuraFlex is intended for use as an alternative to a facemask for achieving and maintaining control of the airway during routine and emergency anaesthetic procedures in fasted patients.

1.2. Intended users and use environment

Medical professionals trained in airway management. AuraFlex is intended to be used in a hospital setting.

1.3. Intended patient population

Adult and pediatric patients from 10 kg and above evaluated as eligible for a supraglottic airway.

1.4. Contraindications

None known.

1.5. Clinical Benefits

Keeps the upper airway open to allow passage of gases.

1.6. Warnings and Cautions

Before insertion, it is essential that all medical professionals using Ambu AuraFlex are familiar with the warnings, precautions, indications, and contraindications found in *Instruction for use*.

warnings /!\

- The product is intended to be used by medical professionals trained in airway management only.
- Always visually inspect the product and perform a functional test after unpacking and prior to use according to section 3.1 Preparation before use, as defects and foreign matters can lead to no or reduced ventilation, mucosal damage, or infection of the patient. Do not use the product if any steps in Preparation before use fails.
- Do not re-use AuraFlex on another patient as it is a single use device. Re-use of a contaminated product can lead to infection.
- 4. AuraFlex does not protect the trachea or lungs from the risk of aspiration.
- Do not use excessive force when inserting and removing AuraFlex as this can lead to tissue trauma.
- The cuff volume or pressure may change in the presence of nitrous oxide, oxygen, or other medical gases which can lead to tissue trauma.
 Ensure to monitor cuff pressure continuously during the surgical procedure.
- Do not use AuraFlex in the presence of lasers and electrocautery equipment as this could lead to airway fire and tissue burns.
- Do not perform direct intubation through AuraFlex as the endotracheal (ET) tube may get stuck leading to insufficient ventilation.

- 9. Ambu AuraFlex is MR unsafe.
- In general, AuraFlex should only be used in patients who are profoundly unconscious and will not resist insertion.
- 11. The overall complication rate for laryngeal mask is low, but the user must exercise professional judgement when deciding whether the use of a laryngeal mask will be appropriate.

The following patients are at higher risk of serious complications including aspiration and inadequate ventilation:

- Patients with upper airway obstruction.
- Non-fasted patients (including those cases where fasting cannot be confirmed).
- Patients suffering from upper gastrointestinal issues (e.g., esophagectomy, hiatal hernia, gastroesophageal reflux disease, morbid obesity, pregnancy > 10 weeks).
- · Patients requiring high pressure ventilation.
- Patients who present with pharyngeal/laryngeal pathology potentially complicating anatomical fit of the mask (e.g., tumors, radiotherapy to the neck involving the hypopharynx, severe oropharyngeal trauma).
- Patients with inadequate mouth opening to permit insertion.

CAUTIONS

- Do not soak, rinse, or sterilize this device as these procedures may leave harmful residues or cause malfunction of the device. The design and material used are not compatible with conventional cleaning and sterilization procedures.
- Before use, always check for compatibility between AuraFlex and the external device to avoid the use of devices not being able to pass through lumen of AuraFlex.
- The cuff pressure should be kept as low as possible while still providing sufficient seal and should not exceed 60 cmH₂O.
- Any signs of airway problems or inadequate ventilation must be monitored regularly and AuraFlex must be repositioned, reinserted or replaced as required to maintain a patent airway.
- 5. Always reconfirm the patency of the airway after any change in the patient's head or neck position.

1.7. Potential adverse events

The use of laryngeal masks is associated with minor adverse effects (e.g., sore throat, bleeding, dysphonia, dysphagia) and major adverse effects (e.g., regurgitation/aspiration, laryngospasm, nerve injury).

1.8. General notes

If, during the use of this device or as a result of its use, a serious incident has occurred, please report it to the manufacturer and to your national authority.

2.0. Device description

AuraFlex is a sterile, single-use laryngeal mask consisting of a flexible patient tube with an inflatable cuff at the distal end. The cuff can be inflated through the check valve allowing the pilot balloon to indicate the inflation/deflation status. The cuff conforms to the contours of the hypopharynx and with its lumen facing the laryngeal opening of the patient. The tip of the cuff presses against the upper oesophageal sphincter and the proximal end of the cuff rests against the base of the tongue.

AuraFlex comes in 6 different sizes. The main components of AuraFlex are seen in figure ①.

Figure 1 (page 7): Overview of AuraFlex parts: 1. Connector; 2. Patient tube; 3. Cuff; 4. Check valve; 5. Pilot balloon; 6. Pilot tube; 7. Nominal length of internal ventilatory pathway*.

* See Table 1 for the nominal length provided in centimeters.

Figure 2 (page 7): Correct position of AuraFlex in relation to AuraFlex parts and anatomical landmarks

AuraFlex parts: 1. Inflatable cuff; 2. Size marking;
3. Ventilatory opening; 4. Ventilatory pathway;
5. Normal depth of insertion marks; 6. Machine end.

Anatomical landmarks: A. Esophagus; B. Trachea; C. Cricoid ring; D. Thyroid cartilage; E. Vocal cords; F. Laryngeai inlet; G. Epiglottis; H. Hyoid bone; I. Tongue; J. Buccal cavity; K. Nasopharynx; L. Incisors.

COMPATIBILITY WITH OTHER DEVICES/EQUIPMENT

AuraFlex can be used in conjunction with:

- Ventilation equipment; 15 mm conical connectors in compliance with ISO 5356-1.
- · Airway management devices; Bronchoscopes*
- Other accessories; Standard 6 % conical Luer syringe, Manometer with standard 6 % conical Luer connector, Water-based lubrication, Suction catheter.

When using instruments through the mask ensure that the instrument is compatible and well lubricated before insertion.

* See Table 1 for information on maximum instrument size.

3.0. Product use

3.1. Preparation before use SIZE SELECTION

Ambu AuraFlex comes in different sizes for use in patients of different weights.

For pediatric patients, it is recommended that Ambu AuraFlex is used by a medical professional familiar with pediatric anesthesia.

See selection guidelines and max. intracuff pressure in Table 1, section 4.0. (Specifications).

INSPECTION OF AURAFLEX

Always wear gloves during the preparation and insertion of Ambu AuraFlex to minimize contamination. Handle AuraFlex carefully as it can be torn or punctured. Avoid contact with sharp or pointed objects.

Check that the pouch seal is intact before opening and discard Ambu AuraFlex if the pouch seal has been damaged.

Closely examine AuraFlex for any damage, such as perforation, scratches, cuts, tears, loose parts, sharp edges etc.

Make sure that the cuff protector has been removed from the cuff.

Check that the interior of the patient tube and the cuff are free from blockage and any loose parts. Do not use AuraFlex if it is blocked or damaged.

Deflate the cuff of AuraFlex completely. Once deflated, check the cuff thoroughly for any wrinkles or folds.

Inflate the cuff to the volume as specified in Table 1. Check that the inflated cuff is symmetrical and smooth. There should not be any bulge nor any sign of leakage in the cuff, pilot tubing or pilot balloon. Deflate the cuff again before insertion.

3.2. Preparations for Use PRE-INSERTION PREPARATION

- Deflate the cuff completely so that the cuff is flat and free of wrinkles by pressing the cuff down onto a flat sterile surface (e.g., a piece of sterile gauze) while at the same time deflating the device with a syringe. 3
- Lubricate the posterior tip of the cuff prior to insertion by applying a sterile, water-based lubricant to the distal posterior surface of the cuff.
- Always have a spare Ambu AuraFlex ready for use.
- Pre-oxygenate and use standard monitoring procedures.
- Check that the level of anesthesia (or unconsciousness) is adequate before attempting insertion. The insertion should be successful at the same level of anesthesia that would be suitable for tracheal intubation.
- The head of the patient should be positioned extended with flexion of the neck in a position normally used for tracheal intubation (i.e., "the sniffing position").

3.3. Insertion

- · Never use excessive force.
- Hold the Ambu AuraFlex like a pencil, with the index finger placed at the transition between the cuff and the patient tube. Your other hand should be placed under the patient's head.
- Insert the tip of the cuff pressing upwards against the hard palate and flatten the cuff against it. (5)
- Verify that the tip of the cuff is flattened against the palate before proceeding – push the jaw gently downwards with your middle finger to open the mouth further.
- Ensure that the tip of the cuff avoids entering the valleculae or the glottic opening and does not become caught up against the epiglottis or the arytenoids. The cuff should be pressed against the patient's posterior pharyngeal wall.
- When the mask is in place, resistance will be felt.
- Before pulling back the hand holding the mask, the other hand is used to press down on the airway tube. This ensures that the Ambu AuraFlex stays in place when the finger is removed.
- After insertion ensure lips are not trapped between patient tube and teeth to avoid trauma to lips.

INSERTION PROBLEMS

- For pediatric patients, a partial rotational technique is recommended in case of placement difficulties.
- Coughing and breath-holding during Ambu AuraFlex insertion indicates inadequate depth of anesthesia – Immediately deepen anesthesia with inhalational or intravenous agents and initiate manual ventilation.
- If you cannot open the patient's mouth sufficiently to insert the mask, check that the patient is adequately anesthetized. Ask an assistant to pull the jaw downwards thus making it easier to see into the mouth and verify the position of the mask.
- For difficulty in maneuvering the angle at the back
 of the tongue when inserting AuraFlex, press the
 tip against the palate throughout or else the tip
 may fold on itself or meet an irregularity in the
 posterior pharynx, e.g., hypertrophied tonsils.
 Should the cuff fail to flatten or begin to curl over
 as it is inserted, withdraw the mask and reinsert
 it. In case of tonsillar obstruction, a diagonal
 movement of the mask is recommended.

3.4. Fixation

Place a Boyle-Davis gag if needed and fixate the Ambu AuraFlex according to surgical needs. It is recommended to use a gauze bite block.

3.5. Inflation

- Without holding the tube, inflate the cuff with just enough air to obtain a seal, equivalent to intracuff pressures of a maximum of 60 cmH₂O. ⑥ Often only half of the maximum volume is sufficient to achieve a seal – please refer to Table 1 for maximum intracuff volumes.
- Monitor the cuff pressure continuously during the surgical procedure with a cuff pressure gauge. This is especially important during prolonged use or when nitrous oxide gases are used.
- Look for the following signs of correct placement:
 The possible slight outward movement of the tube
 upon cuff inflation, the presence of a smooth oval
 swelling in the neck around the thyroid and cricoid
 area, or no cuff visible in the oral cavity.
- The mask may leak slightly for the first three or four breaths before settling into position in the pharynx. In case leakage persists, check that there is adequate depth of anesthesia and that the pulmonary inflation pressures are low before assuming that reinsertion of AuraFlex is necessary.

3.6. Verification of correct position

 Correct placement should produce a leak-free seal against the glottis with the tip of the cuff at the upper oesophageal sphincter.

- The vertical line on the patient tube should be oriented anteriorly towards the patient's nose.
- AuraFlex is inserted correctly when the patient's incisors are between the two horizontal lines on the patient tube. (2), item 5. Reposition the mask if the patient's incisors are outside this range.
- The position of AuraFlex can be assessed by capnography, by observation of changes in tidal volume (e.g., a reduction in expired tidal volume), by auscultating bilateral breath sounds and an absence of sounds over the epigastrium and/or by observing chest rise with ventilation. If you suspect that AuraFlex has been positioned incorrectly, remove and reinsert – and ensure that anesthetic depth is adequate.
- Visual confirmation of anatomically correct position is recommended, e.g., by using a flexible scope.

UNEXPECTED REGURGITATION:

- Regurgitation may be caused by inadequate level of anesthesia. The first signs of regurgitation may be spontaneous breathing, coughing or breath-holding.
- If regurgitation occurs, if oxygen saturation remains at acceptable levels, AuraFlex should not be removed. This should be managed by putting the patient in a "head-down" position. Briefly

disconnect the anesthetic circuit so that the gastric contents are not forced into the lungs. Check that anesthetic depth is adequate and deepen anesthesia intravenously, if appropriate.

 Apply suction through the mask's patient tube and through the mouth. Suction the tracheobronchial tree and inspect the bronchi using a flexible scope.

3.7. Usage with other devices/equipment ANESTHETIC SYSTEM AND VENTILATION BAG

The mask can be used for either spontaneous or controlled ventilation.

During an esthesia, nitrous oxide may diffuse into the cuff causing an increase in cuff volume/pressure. Adjust cuff pressure just enough to obtain an adequate seal (cuff pressure should not exceed 60 cm H_2O).

The anesthetic breathing system must be adequately supported when connected to AuraFlex to avoid rotation of the mask.

USAGE WITH SPONTANEOUS VENTILATION

AuraFlex is suitable for spontaneously breathing patients when used with volatile agents or intravenous anesthesia on condition that anesthesia is adequate to match the level of surgical stimulus and the cuff is not overinflated.

USAGE WITH POSITIVE PRESSURE VENTILATION

When applying positive pressure ventilation, ensure that the seal is adequate. To improve the seal the following is suggested:

- Optimize placement of AuraFlex by head turning or traction.
- Adjust the cuff pressure. Try both lower and higher pressures (a poor cuff seal may be caused by either too low or too high cuff pressure).
- If leakage should occur around the cuff, remove the mask and reinsert while ensuring that anesthetic depth is adequate.

MAGNETIC RESONANCE IMAGING (MR)

AuraFlex is MR-unsafe. The device presents a projectile hazard.

3.8. Removal procedure

Removal should always be carried out in an area where suction equipment and the facility for rapid tracheal intubation are available.

Do not remove AuraFlex with the cuff fully inflated to prevent tissue trauma and laryngospasm.

3.9. Disposal

Dispose used Ambu AuraFlex in a safe manner according to local procedures.

4.0. Specifications

Ambu AuraFlex is in conformity with ISO 11712 Anaesthetic and respiratory equipment – Supralaryngeal airways and connectors.

	Ped	iatric	Adult					
Mask Size	#2	#21/2	#3	#4	#5	#6		
Patient weight	10 – 20 kg	20 – 30 kg	30 – 50 kg	50 – 70 kg	70 – 100 kg	> 100 kg		
Maximum intracuff volume	10 ml	14 ml	20 ml	30 ml	40 ml	50 ml		
Maximum intracuff pressure	60 cmH₂O							
Airway connector	15 mm male (ISO 5356-1)							
Maximum instrument size*	5.5 mm	5.5 mm	7.0 mm	7.0 mm	8.0 mm	8.0 mm		
Inflation Valve Luer cone compatibility	Luer cone compatible with ISO 594-1 and ISO 80369-7 compliant equipment							
Appropriate storage condition	10 °C (50 °F) to 25 °C (77 °F)							
Approximate mask weight	22 g	25 g	34 g	39 g	49 g	54 g		
Internal volume of ventilatory pathway	8.0 ± 0.3 ml	8.6 ± 0.8 ml	13.2 ± 0.8 ml	13.6 ± 1.2 ml	17.6 ± 1.5 ml	18.5 ± 0.8 ml		
Pressure drop as determined according to ISO 11712 annex C	0.8 cmH₂O at 30 l/min	0.7 cmH₂O at 30 l/min	0.8 cmH₂O at 60 l/min	0.7 cmH₂O at 60 l/min	0.6 cmH₂O at 60 l/min	0.6 cmH₂O at 60 l/min		
Min. Interdental gap	17 mm	20 mm	21 mm	24 mm	27 mm	29 mm		
Nominal length of the Internal ventilatory pathway	21.2 ± 1.3 cm	21.2 ± 1.3 cm	23.7 ± 1.4 cm	23.7 ± 1.4 cm	26.2 ± 1.6 cm	26.2 ± 1.6 cm		

Table 1: Specifications for Ambu AuraFlex.

A full list of symbol explanations can be found on https://www.ambu.com/symbol-explanation

No part of this documentation may be reproduced in any form, including photocopying, without the prior written permission of the copyright owner.

^{*}The maximum instrument size is intended as a guide for selecting the appropriate diameter of a device to be passed through the patient tube of AuraFlex.

[©] Copyright 2021 Ambu A/S, Denmark. All rights reserved.

Ambu



Ambu A/S

Baltorpbakken 13 2750 Ballerup Denmark T +45 72 25 20 00

ambu.com

UK RP

0086

Ambu Ltd

First Floor, Incubator 2

Alconbury Weald Enterprise Campus

Alconbury Weald

Huntingdon PE28 4XA

United Kingdom

www.ambu.co.uk

