

Ambu® Perfit ACE® & Mini Perfit ACE®



Key Benefits

- Stores completely flat
- Radio translucent
- Sizing line for easy and accurate sizing
- Standard size markings on both Adult and Mini
- Chin flip feature to ease airway management
- Nasal cannula tubing holder support
- Facilitates intubation
- Latex free
- Disposable
- When shipped, the default position of the the Adult Ace collar is in the Neckless position, and the default position of the Mini Ace is the Infant position

Ambu® Perfit ACE® & Mini Perfit ACE®

28 sizes in two collars

The Ambu Perfit ACE is a unique adjustable collar. Combined with the Mini Perfit ACE these two collars offer an incredible 28 sizes ensuring that patients can now have the most effective and accurate fit ever.

Made from radio translucent materials, the collar is packaged completely flat and has a unique chin-flip feature to ease airway management. Ambu Perfit ACE is a disposable adult adjustable collar for immobilizing the cervical column.

Ambu Perfit Ace is an adjustable collar designed to fit any adult with 16 different settings ranging from Neckless to Tall sizes.

The Ambu Perfit ACE incorporates latches that allow the rescuer to size the collar exactly to the victim's neck size. Two safety lock buttons hold the collar in place when adjusted.

Ambu Mini Perfit ACE is an adjustable collar designed to fit children and small adults, with 12 different settings to provide an accurate fit.

Specifications


Ambu Perfit ACE	
Adjustable sizes	Neckless, Short, Regular and Tall
Number of settings	16 adjustments
Dimensions	554 x 215 x 14 mm (21.8" x 8.5" x 0.6")
Weight	200 g (7 ounces)
Color	White/light grey
Ambu Mini Perfit ACE	
Adjustable sizes	Infant, Pediatric and Small Adult
Number of settings	12 adjustments
Dimensions	447 x 188 x 14 mm (17.6" x 7.4" x 0.6")
Weight	100 g (3.5 ounces)
Color	Yellow/light grey




Ambu® Perfit ACE®



Ambu® Mini Perfit ACE®

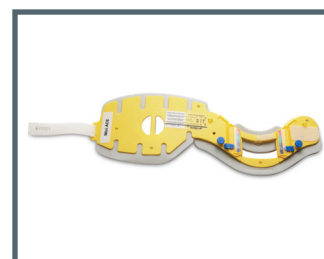

Ambu® Perfit ACE®
Ambu® Mini Perfit ACE®

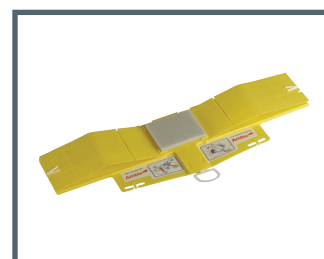

Ambu® Head Wedge®

Material

Product	Material
Ambu Perfit ACE	
White/light grey plastic	Polypropylene
Foam	Polyethylene
Ambu Mini Perfit Ace	
Yellow plastic	Polypropylene
Foam	Polyethylene

Ordering specifications

Item no.	Description
000 281 000	Perfit ACE Extrication Collar, 30/Case
000 281 151	Perfit ACE packed with Head Wedge CID, 30/Case
000 281 106	Mini Perfit ACE Extrication Collar, 30/Case
000 281 161	Mini Perfit ACE packaged with Head Wedge CID, 30/Case

Available packaged with the Ambu® Head Wedge®.

Ambu Australia Pty Ltd
Customer Services
 Tel.: 1300 233 118
 Fax. 1300 156 526
 aussales@ambu.com
 www.ambuaustralia.com.au

Ambu Asia (Sales office)
Customer Services
 Tel.: +604 252 9100
 Fax. +604 644 8932
 csasia@ambu.com