

Ambu® Redi-ACE Adult & Redi-ACE Mini



Key Benefits

- **Easy to use**

Ambu Redi-ACE has an automatic flip up chin and two locking buttons. It can be placed on the patient using only 3 steps; sizing – locking – placement.

- **2 collars for all patients**

With 8 sizes of Ambu Redi-ACE and 6 sizes of Ambu Redi-ACE Mini the majority of patients can be treated.

- **Practical characteristics**

Ambu Redi-ACE can be stored completely flat. Additionally the Ambu Redi-ACE has cannula hooks to secure the oxygen tube.

Ambu® Redi-ACE Adult & Redi-ACE Mini

The Ambu Redi-ACE is a unique adjustable collar. Combined with the Redi-ACE Mini these two collars offer 14 sizes ensuring that patients can have the most effective and accurate fit possible.

Made from radio translucent materials, the collar is packaged completely flat and has a unique automatic flip up chin.

Ambu Redi-ACE is a disposable adjustable collar for immobilizing the cervical column.

The Ambu Redi-ACE incorporates latches that allow the rescuer to size the collar as accurate as possible to the victim's neck size. Two safety lock buttons hold the collar in place when adjusted.

Specifications

Ambu Redi-ACE	
Adjustable sizes	Neckless, Short, Regular and Tall
Number of settings	8 adjustments
Dimensions	554 x 215 x 14 mm (21.8" x 8.5" x 0.6")
Weight	200 g (7 ounces)
Color	Blue/light grey
Translucency	CT and MRI compatible
Storage temperature	-30° to +70° C
Operation temperature	-10° to +40° C
Ambu Redi-ACE Mini	
Adjustable sizes	Infant, Pediatric and Small Adult
Number of settings	6 adjustments
Dimensions	447 x 188 x 14 mm (17.6" x 7.4" x 0.6")
Weight	100 g (3.5 ounces)
Color	Yellow/light grey
Translucency	CT and MRI compatible
Storage temperature	-30° to +70° C
Operation temperature	-10° to +40° C




Ambu Redi-ACE Adult



Ambu Redi-ACE Mini

Material

Product	Material
Ambu Redi-ACE	
Blue/light grey plastic	Polypropylene
Foam	IXPE
Ambu Redi-ACE Mini	
Yellow/light grey plastic	Polypropylene
Foam	IXPE

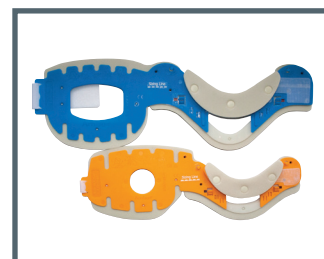

Flat storage

Ordering specifications

Item no.	Description
472 001 000	Redi-ACE Extrication Collar, 30/Case
472 002 000	Redi-ACE Mini Extrication Collar, 30/Case



Locking buttons



Cannula hooks

Ambu Australia Pty Ltd
Customer Services
 Tel.: 1300 233 118
 Fax. 1300 156 526
 aussales@ambu.com
 www.ambuaustralia.com.au

Ambu Asia (Sales office)
Customer Services
 Tel.: +604 252 9100
 Fax. +604 644 8932
 csasia@ambu.com