

## DATASHEET

AmbuMan® Basic is designed to fulfill all requirements for the BLS training algorithm.

With the new iQF App, you can provide training to a high level of precision.


### **CPR becomes the most important factor when a person suffers from cardiac arrest.**

After three minutes without oxygen, the brain cells start to die.<sup>1</sup>

### **iQF, Training that Works for Life.**

Helping instructors to train. The iQF application provides data-driven information and accurate measurements of relevant CPR parameters.

Helping to ensure students' CPR technique is correct in a fun, informative and interactive way.

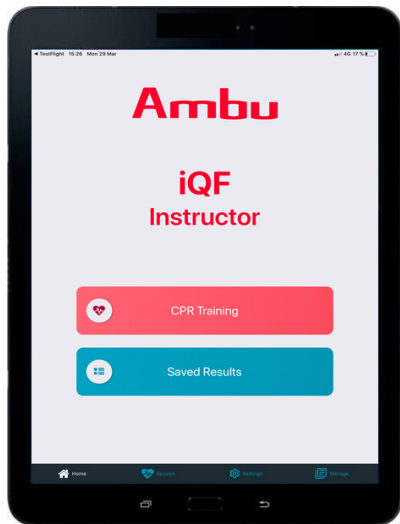

Ambu iQF Instructor interface

### **Realistic Ventilation and Compression**

The AmbuMan Basic lets students see, hear and feel the ventilation.

# AmbuMan® Basic iQF



## KEY BENEFITS

- Robust, reliable and cost efficient BLS training manikin with anatomically correct closed torso
- Inbuilt visual indicator, providing instant feedback of compression amplitude
- Ambu hygienic system for mouth-to-mouth and mouth-to-nose ventilation to help prevent cross-contamination
- Realistic ventilation
- Enables easy AED Training

## TRAIN IQF, TRAIN FOR LIFE

Real time visualisation of the compression depth and the correct ventilation volume at your fingertips with the Instructor and Trainee iQF App. Connected by Bluetooth®.

# Ambu

## SPECIFICATIONS

|                           |         |
|---------------------------|---------|
| Weight Torso (single)     | 6 kg    |
| Length                    | 70 cm   |
| Width                     | 37 cm   |
| Weight 4 x Torsos and Bag | 30 kg   |
| Max. ventilation volume   | 1200 ml |
| Max. compression depth    | 7 cm    |

## MATERIALS

|                    |                                       |
|--------------------|---------------------------------------|
| Basic unit         | ABS (Acrylonitrile Butadiene Styrene) |
| Head               | Polyvinyl chloride, hard              |
| Chest skin         | Polyvinyl chloride, soft              |
| Face piece         | Polyvinyl chloride, soft              |
| Headbag disposable | Polyethylene                          |
| Training suit      | Polyester                             |
| Carry bag          | PVC Coated Nylon                      |

## ORDERING SPECIFICATIONS

| Item no.   |   |
|------------|---|
| A268008000 | AmbuMan Basic iQF incl. 1 AmbuMan Basic Torso with training jacket, carry case, 1 face piece and 25 head bags               |
| A268009000 | AmbuMan Basic iQF 4 Pack iQF incl. 4 AmbuMan Basic Torsos with training jacket, carry case, 4 face pieces and 100 head bags |

| Accessories |   |
|-------------|---|
| 234 000 702 | Ambu head bags, disposable (set of 100) |
| 234 000 703 | Ambu face piece (set of 5)              |
| 234 000 725 | Ambu Guedel face piece                  |


\* Bluetooth is a registered trademark of Bluetooth SIG, inc.



Compression Indicator



Ambu's unique Hygienic System

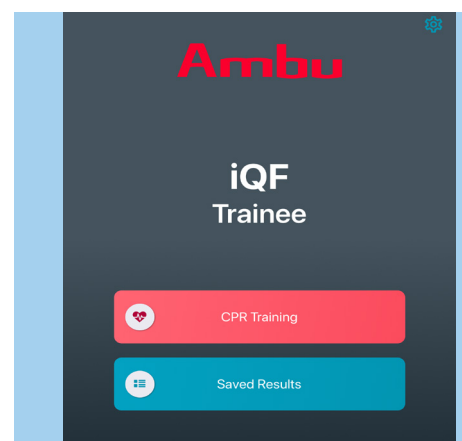

Ambu iQF Trainee interface

# Ambu

### Ambu Australia Pty Ltd

Unit 5, 2 Daydream Street, Warriewood, NSW 2102  
Tel: 1300 233 118 | Email: aussales@ambu.com  
www.ambuaustralia.com.au  
AA22 0016 V01

### References

1. Intensive Care NSW, Cardiac Arrest, viewed 20 June 2022, <<https://aci.health.nsw.gov.au/networks/icnsw/patients-and-families/patient-conditions/cardiac-arrest>>.